



Image from jax.org

# Mouse Strain Datasheet

## Strain Name

Pkd1ko

**MGI Gene ID**

MGI:3617325

**Full Allele Name**

Pkd1<tm1Ggg>

**Type of Allele**

knockout allele

**Human Gene (HGNC)**

PKD1

**Genetic Background**

C57BL/6J

**Commercial Source**

n/a

**Stock Number**

n/a

**Link**

n/a

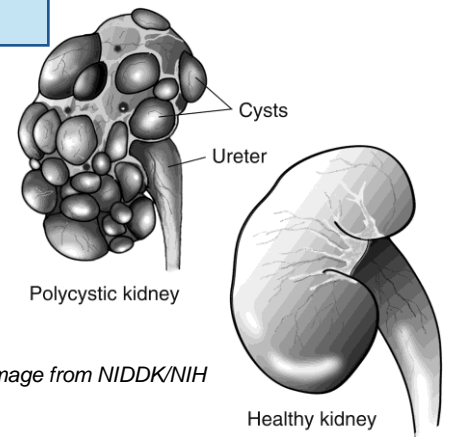


Image from NIDDK/NIH

**Genotyping Protocol**

<https://pubmed.ncbi.nlm.nih.gov/15579506/>

### Strain Details

These mice possess a true null allele and lacked expression of the targeted gene. Mice that are homozygous for this allele are embryonically lethal.

### Validation or publication

<https://pubmed.ncbi.nlm.nih.gov/15579506/>

**Contact Name**

Patricia Outeda

**Email**

pgarcia@som.umaryland.edu

## PCR Protocol for Genotyping:

Pkd1ko

### A. Digestion of mouse tail or ear notch, and **embryo tail (in red)**:

1. Add 100µL of Tissue Digestion Buffer and 2µL of Proteinase K per tail (~1-2mm length). For embryos tail add **50µL of Tissue Digestion Buffer and 1µL of Proteinase K**. Make sure tail is immersed in the buffer.
2. In a thermocycler incubate at 55°C for 1 h followed by 95°C for 8 min to inactivate the enzyme and hold at 10°C. For embryos incubate at 55°C for 30 min followed by 95°C for 8 min and hold at 10°C.
3. Vortex and store at 4°C (-20°C for long storage) or use immediately to set up the PCR.

### B. PCR Genotyping Protocol

Primers			
KO forward 1	5'-	AAT AGG GGT GGG GCT TGT GGG TCG	-3'
KOK back 1	5'-	TGG CGA AAG GGG GAT GTCG CTG C	-3'
M3.2B	5'-	TAC TCA CAC CTC CAC CAG TGC	-3'
	5'-		-3'

PCR Reaction		PCR Conditions		
BioMix (Bioline)	10.0 µL		Heated Lid	105°C
Primers (@10 µM each)	0.8 µL		Initial Denaturation	94°C 5 min
			Number of Cycles	x35
ddH <sub>2</sub> O	7.2 µL		94°C	20 sec
			59°C	35 sec
			72°C	35 sec
DNA template	2.0 µL		Final Extension	72°C 10 min
Total Volume	20.0 µL		Final Hold	10°C

PCR Product Size (bp)	
Wild type band	180 bp
ko band	250bp

### C. Reagents

Reagent	Cat #	Final Concentration	Working Concentration
<b>Tissue Digestion Buffer for ear notch or tail</b>			
Tris pH8.5		50mM	
EDTA		1mM	
Tween20		0.5%	
<b>Proteinase K (Invitrogen)</b>	25530-015	20mg/mL	
<b>BioMix (Bioline)</b>	BIO-25012		